Concept House Cardiff Road Newport South Wales NP10 8QQ

ne undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) he Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of Comptroller-General, hereby certify that annexed hereto is a true copy of the documents riginally filed in connection with patent application GB0301278.8 filed on 21 January 3.

cordance with the Patents (Companies Re-registration) Rules 1982, if a company named s certificate and any accompanying documents has re-registered under the Companies 980 with the same name as that with which it was registered immediately before reration save for the substitution as, or inclusion as, the last part of the name of the words c limited company" or their equivalents in Welsh, references to the name of the my in this certificate and any accompanying documents shall be treated as references to me with which it is so re-registered.

rdance with the rules, the words "public limited company" may be replaced by p.l.c., c. C. or PLC.

stration under the Companies Act does not constitute a new legal entity but merely the company to certain additional company law rules.

Signed

Dated 11 May 2007

dti ADTI SERVICE

atents Form 1/77 Patents Act 1977 (Rule 16) The Patent Office Request for grant of a patent 21 JAN 2003 (See the notes on the back of this form. You can also get an Cardiff Road explanatory leaflet from the Patent Office to help you fill in Newport this form) **NEWPORT** South Wales NP10 8QQ Your Reference P51255D GB Patent Application Number 0301278.8 21 JAN 2003 (The Patent Office will fill in this part) Fill name, address and postcode of the or of Elekta AB (publ) each applicant (underline all surnames) P O Box 7593 S-103 93 Stockholm 851582700 Patents ADP number (if you know it) If the applicant is a corporate body, give the Sweden country/state of its incorporation 4. Title of the invention Computed Tomography Scanning Name of you agent (if you have one) Fry Heath & Spence LLP The Gables "Address for service" in the United Kingdom Massetts Road to which all correspondence should be sent (including the postcode) Horley Surrey RH6 7DQ Patents ADP number (if you know it) 05880273001 84 20224001 If you are declaring priority from one or more Country Priority application number Date of filing earlier patent applications, give the country (if you know it) (day / month / year) and the date of filing of the or each of these earlier applications and (if you know it) the or each application number. If this application is divided or otherwise Number of earlier application Date of filing derived from an earlier UK application, (day / month / year) give the number and the filing date of the earlier application. Is a statement of inventorship and of right Yes to grant of a patent required in support of

this request? (Answer 'Yes' if:

applicant, or

See note (d))

any applicant named in part 3 is not an inventor, or there is an inventor who is not named as an

any named applicant is a corporate body.

Patents Form 1/77


Enter the number of sheets for any of the

following items you are filing with this form. Do not count copies of the same document Continuation sheets of this form Nil Description Claims(s) Abstract Drawing(s) 10. If you are also filing any of the following, state how many against each item. Priority documents Translations of priority documents Statement of inventorship and right Two to grant of a patent (Patents Form 7/77) Request for preliminary examination One and search (Patents Form 9/77) Request for substantive examination (Patents Form 10/77) Any other documents (please specify) 11. I/We request grant of a patent on the basis of this application. Signature Date 20 January 2003 12. Name and daytime telephone number of Michael Downing person to contact in the United Kingdom 01865 841060 Warning After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you filing or applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked. Notes If you need help to fill in this form or you have any questions, please contact the Patent Office on 0645 500505. a) b) Write your answers in capital letters using black ink or you may type them. If there is not enough space for all the relevant details on any part of this form, please continue on a separate c) sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form. d) If you have answered 'Yes' Patents Form 7/77 will need to be filed. Once you have filled in the form you must remember to sign and date it. For details of the fee and ways to pay please contact the Patent Office.

Patents Form 7/77

Patents Act 1977 (Rule 15)

Statement of inventorship and of right to grant of a patent

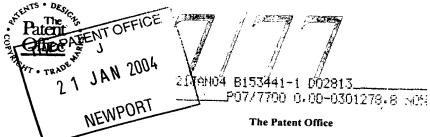

Cardiff Road Newport South Wales NP10 8QQ

1.	Your Reference	P51255D GB
2.	Patent Application Number (if you know it)	0301278.8
3.	Fill name of the or of each applicant	Elekta AB (publ) P O Box 7593 S-103 93 Stockholm Sweden
4.	Title of the invention	Computed Tomography Scanning
5.	State how the applicant(s) derived the right from the inventor(s) to be granted a patent	BROWN: By virtue of employment JAFFRAY & SIEWERDSEN: by virtue of a research agreement dated 5 April 2002 VAN HERK & SONKE: by virtue of a research agreement dated 27 September 1999
5.	How many, if any, additional Patents Forms 7/77 are attached to this form? (see note (c))	One additional page
·.		I/We believe that the person(s) named over the page (and on any extra copies of this form) is/are the inventor(s) of the invention which the above patent application relates to. Signature Date 19 January 2004
	Name and daytime telephone number of person to contact in the United Kingdom	Michael Downing 01865 841060

Notes

- a) If you need help to fill in this form or you have any questions, please contact the Patent Office on 0645 500505.
- b) Write your answers in capital letters using black ink or you may type them.
- c) If there are more than three inventors, please write the names and addresses of the other inventors on the back of another Patents Form 7/77 and attach it to this form.
- d) When an application does not declare any priority, or declares priority from an earlier UK application, you must provide enough copies of this form so that the Patent Office can send one to each inventor who is not an applicant.
- e) Once you have filled in the form you must remember to sign and date it.

Enter the full names, addresses and postcodes of the inventors in the boxes and underline the surnames

Kevin John BROWN

200 New Street

Horsham

West Sussex

RH13 5EJ

760603600

Patents ADP number (if you know it)

David JAFFRAY

81 Poplar Heights Drive

Etobicoke

Ontario

Canada

M9A4Z3

8796195001

Patents ADP number (if you know it)

Jeffrey SIEWERDSEN

230 Queen's Quay W. #1215

Toronto

Ontario

Canada

M5J2Y7

Reminder Have you signed the form? 979670300

Patents ADP number (if you know it)

Patents Form 7/77

Enter the full names, addresses and postcodes of the inventors in the boxes and underline the surnames

Marcel VAN HERK
JF Berghoefplantsoen 11
1064 DE Amsterdam
The Netherlands

8796257001

Patents ADP number (if you know it)

Jan-Jakob SONKE Lijnbaansgracht 118 I 1016 VV Amsterdam The Netherlands

9797268001

Patents ADP number (if you know it)

Reminder Have you signed the form?

Patents ADP number (if you know it)

Patents Form 7/77

Computed Tomography Scanning

FIELD OF THE INVENTION

The present invention relates to scanning by computed tomography (CT).

BACKGROUND ART

CT scanning is a process for imaging the internal structure of a patient. In conventional CT scanning, a beam of x-rays is projected through the patient and its attenuation is measured. At the same time, the apparatus is rotated about an axis passing longitudinally through the patient. Thus, data is acquired as to the attenuation of the beam in each direction in the plane in which rotation takes place. From this data, the internal structure of the patient on that plane can be computed. The patient or apparatus is then indexed along the axis and a further plane (known as a 'slice') is then investigated. A three dimensional image of the patient can then be constructed from the various slices.

One problem is that over the time required to acquire the necessary slices, the patient is not motionless. Gross motor movement can be avoided by suitable instruction to the patient, but even so each slice is acquired at a

different phase of the breathing cycle. This results in a beating artefact due to the different frequency of breathing and slice acquisition.

Two ways have been used to solve this problem. On is to trigger the CT on a particular phase of the patients breathing. This is termed 'respiration gated CT' and implies that one CT slice is acquired for every breath. This means that it takes a long time to acquire a complete volume of data.

Another technique is to monitor the phase of the patients breathing whilst acquiring CT slices continuously. Once the data is acquired, slices that have comparable breathing phase are selected from the complete set and these are then used to visualise the volume. This has the advantage that any phase can be selected retrospectively and therefore the effect of breathing can be studied. This is termed 'respiration correlated CT'.

Conventional CT scans have the disadvantage that the resolution along the axis is poor since it corresponds to the slice thickness. It is theoretically straightforward to increase this, but doing so results in a correspondingly longer acquisition time, or the need to rotate the apparatus correspondingly faster. Both options also give rise to an attendant reduction in contrast in the measured beam. Accordingly, 'cone beam CT' methods have been developed, in which a conical beam of radiation is directed at the patient and a two-dimensional image acquired via a flat panel detector. This apparatus is then rotated around the patient axis and a three-dimensional image is reconstructed from the set of two-dimensional images. As the individual slices are eliminated, there is the same resolution in all directions of the image. Likewise, as there are no slices the above-mentioned breathing artefact is absent since there can be no variation in patient position between slices.

SUMMARY OF THE INVENTION

We have however found that there are other artefacts in the reconstructed volume data of cone beam CT systems, which we have traced to patient breathing movements. In addition, the motion is not measurable in the

reconstructed volume data. This can be a particular problem in cone beam systems due to the long time required for acquisition, typically 1-2 minutes.

The techniques used in conventional CT scanners cannot be used directly in a cone beam system as the data is acquired in 2D projection images, and therefore slices cannot be selected from the resulting data. However, respiration correlation techniques could be applied to the acquired projection images rather than the reconstructed CT volume. To achieve this, we propose monitoring the phase of the patients breathing while acquiring projection images continuously. On completion of the acquisition, projection images that have comparable breathing phases can be selected from the complete set, and these are used to reconstruct the volume data using similar techniques to those of conventional CT. An advantage is that any phase can be selected and therefore the effect of breathing can be studied.

Breath control systems are available, intended for use in conventional CT scanning, and which could be used to monitor the patient's breathing. As an alternative, however, it is possible to use a feature in the projection image(s) to determine the breathing phase. A suitable feature is the position of the patient's diaphragm. This can then be used to select the relevant images to be used in the projection process.

It is known in the field of conventional CT scanning to be advantageous to prompt the patient visually and audibly in order to ensure a regular amplitude and pattern of breathing. Techniques such as these could usefully be applied in the present invention.

Furthermore this feature in the projection images can be used to control delivery of therapeutic radiation dependent on the patient's breathing cycle, to ensure that the tumour is in the correct position when the radiation is delivered. This will provide a direct measure of the patient's breathing phase, a significant improvement as compared to current methods that use external markers affixed to the patient. The use of a 3D volume data set generated using the same

patient position and contemporaneous with the treatment will remove significant uncertainties.

The present invention further provides a radiotherapy device comprising a respiration correlated cone beam CT scanner and a source of therapeutic radiation, in which therapeutic radiation is delivered during the scan at times correlated with the patient's breathing cycle.

BRIEF DESCRIPTION OF THE DRAWINGS

An embodiment of the present invention will now be described by way of example, with reference to the accompanying figures in which;

Figure 1 is a view of a cone beam CT scanner according to the present invention, viewed along the axis of rotation thereof;

Figure 2 is a schematic view of the system incorporating such a scanner; and

Figure 3 shows a treatment apparatus including the scanner of the present invention.

DETAILED DESCRIPTION OF THE EMBODIMENT

Figure 1 shows a cone beam CT scanner. A patient 10 is supported on a couch 12 which may be of any suitable design. Couches typically allow the elevation and longitudinal position of the patient to be adjusted and this may be provided for as desired.

An x-ray source 14 is arranged to project a wide beam 16 of radiation generally directed towards the isocentre 18 of the patient. The source 14 is rotatable around the isocentre 18 on a rotational support 20. The support can, for example, be in the form of a ring or annulus around the patient 10 and couch 12 in which the source is mounted, or it can be a C-arm, or any suitable support allowing the source to rotate, or any combination thereof.

A two-dimensional flat-panel detector 22 is also mounted on the support 20, opposite the source 14 and arranged to rotate in synchronism therewith. If the support includes a C-arm then this can be achieved by mounting the detector on the opposite arm.

Thus, radiation emitted by the source 14 is partially absorbed by the patient and the attenuated signal is detected by the flat panel detector 22. The source 14 and detector 22 are then indexed rotationally and a fresh image obtained. This is repeated until sufficient images are acquired to reconstruct the volume data, typically one complete rotation.

Figure 2 shows the system as a whole. The scanner of figure 1 is shown, together with cables linking the source 14, detector 22 and rotational support 20 to a plurality of computing means 24, 26 which process the data generated including the images, source intensity (etc), and rotational support position. Data is output via any suitable means, depicted generally as a monitor 28 but not limited thereto, and the system is controlled by any suitable input means, again depicted generally as a keyboard 30 but likewise not especially limited thereto.

As mentioned above, we have found that there are artefacts in the reconstructed volume data of cone beam CT systems, which we have traced to patient breathing movements. To overcome or alleviate these, respiration correlation techniques are applied to the acquired projection images by the computing means 24, 26. This differs from conventional respiration-correlated CT scanning in acting on the acquired projection images rather than the reconstructed CT volume.

To assist in this process, a breath control system is provided at 32 to monitor the phase of the patients breathing while the projection images are acquired. On completion of the acquisition, projection images that have comparable breathing phases can be selected from the complete set, and these are used to reconstruct the volume data using cone beam CT techniques. As a

result, any phase or range of phases can be selected and therefore the effect of breathing can be studied if desired.

As an alternative to the breath control system, it is possible to use a feature in the projection image(s) to determine the breathing phase, such as the position of the patient's diaphragm. This can then be used to select the relevant images to be used in the projection process.

An alert system including a light 34 and a buzzer 36 is provided, to prompt the patient visually and audibly in order to ensure a regular amplitude and pattern of breathing. Other alerts could of course be employed, such as other forms of visible prompts including (for example) movable devices, and other forms of audible prompts including (for example) speakers, percussive devices or any other form of controllable sound generation apparatus.

Figure 3 shows a system including a therapeutic source of radiation 38 arranged to emit a suitably collimated beam of therapeutic radiation 40. This allows simultaneous scanning and treatment. If the radiation from source 14 continues during the treatment the selected feature (above) in the projection images can be used to control delivery of therapeutic radiation from the source 38, dependent on the patient's breathing cycle. This ensures that the tumour is in the correct position when the radiation is delivered. This will provide a direct measure of the patient's breathing phase, a significant improvement as compared to current methods that use external markers affixed to the patient. The use of the same direct measure of the patient's breathing phase and patient position to generate the 3D volume data set and control the treatment delivery will remove significant uncertainties.

It will of course be understood that many variations may be made to the above-described embodiment without departing from the scope of the present invention.

CLAIMS

- 1. A respiration correlated cone beam CT scanner.
- 2. A method of cone beam CT scanning including respiration correlation.
- 3. A method of cone beam CT scanning in which respiration correlation techniques are applied to the acquired two-dimensional projection images.
- 4. A method of cone beam CT scanning according to claim 3 in which the phase of the patient's breathing is monitored continuously during acquisition of projection images.
- 5. A method of cone beam CT scanning according to claim 4 in which projection images that have comparable breathing phases are selected from the complete data set on completion of the acquisition and are used to reconstruct the volume data.
- 6. A method of cone beam CT scanning according to claim 4 or claim 5 in which a feature in the projection image(s) is used to determine the breathing phase.
- 7. A method of cone beam CT scanning according to claim 6 in which the feature is the position of the patient's diaphragm.
- 8. A method of cone beam CT scanning according to any one of claims 2 to 7 in which visual and/or audible prompts are provided for the patient's breathing.
- 9. A cone beam CT scanner according to claim 1 or a method of cone beam CT scanning according to any one of claims 2 to 8 in which therapeutic radiation is delivered during the scan at times correlated with the patient's breathing cycle.
- A cone beam CT scanner including means for acquiring information as to the patient respiration cycle and means for selection of acquired two-

dimensional projection images from the set of data acquired during a scan on the basis of the respiration cycle information.

- 11. A cone beam CT scanner according to claim 10 adapted to monitor the phase of the patient's breathing continuously during acquisition of projection images.
- 12. A cone beam CT scanner according to claim 11 arranged to select projection images that have comparable breathing phases from the complete data set on completion of the acquisition and to use these to reconstruct the volume data.
- 13. A cone beam CT scanner according to claim 11 or claim 12 including means for detecting a respiration-cycle-correlated feature in the projection image(s) thereby to determine the breathing phase.
- 14. A cone beam CT scanner according to claim 13 in which the feature is the position of the patient's diaphragm.
- 15. A cone beam CT scanner according to any one of claims 9 to 14 including means to provide visual and/or audible prompts for the patient's breathing.
- 16. A radiotherapy device comprising a respiration correlated cone beam CT scanner and a source of therapeutic radiation, in which therapeutic radiation is delivered during the scan at times correlated with the patient's breathing cycle.
- 17. A cone beam CT scanner substantially as described herein with reference to and/or as illustrated in the accompanying figures.
- 18. A method of cone beam CT scanning substantially as described herein with reference to and/or as illustrated in the accompanying figures.

ABSTRACT

COMPUTED TOMOGRAPHY SCANNING

Artefacts in the reconstructed volume data of cone beam CT systems can be removed by the application of respiration correlation techniques to the acquired projection images. To achieve this, the phase of the patients breathing is monitored while acquiring projection images continuously. On completion of the acquisition, projection images that have comparable breathing phases can be selected from the complete set, and these are used to reconstruct the volume data using similar techniques to those of conventional CT. Any phase can be selected and therefore the effect of breathing can be studied. It is also possible to use a feature in the projection image(s) such as the patient's diaphragm to determine the breathing phase. This feature in the projection images can be used to control delivery of therapeutic radiation dependent on the patient's breathing cycle, to ensure that the tumour is in the correct position when the radiation is delivered.

(Fig 1)

